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Dear users,


we congratulate on the buying of a technically matured and the best available technology product.

This instructions helps you and the carers as a basis for sure, safe use of the toilet wheel chair and a sensitive use of the product's possibilities.

The instructions describes the preparation, service, use, servicing and cleaning of the product.

This instructions is to be perused carefully and must be always accessible to all users of the product.

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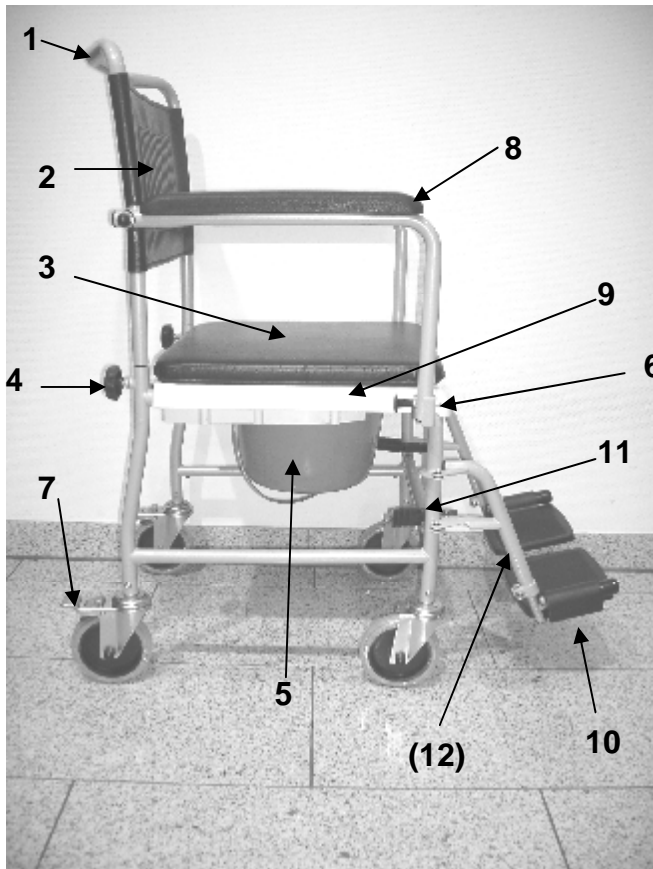Toilet wheel chairs are necessary when using or achievement of the toilet is not possible because of a handicap. Toilet wheel chairs are also usable as room wheel chairs. Toilet wheel chairs can be pushed over customary toilet bowls.

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1. Wheeled Drop-Arm Commode
2. Bucket with cover
3. Footrest
4. Heel tape
5. Operating instructions

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This product is suitable for the reapplication.


With the processing the standards of the manufacturer must be considered. This standards are made available on inquiry.

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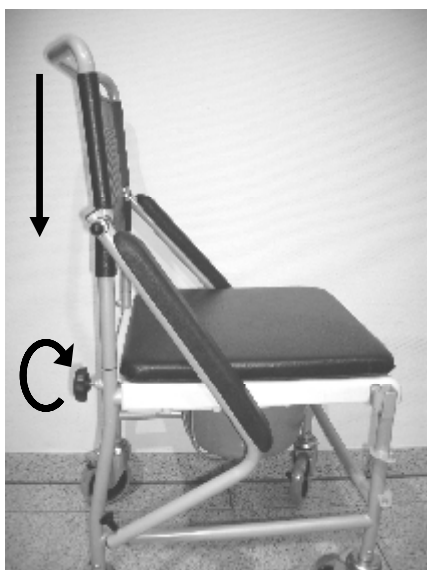The Wheeled Drop-Arm Commode is mountable nearly without tools

Put the premounted backrest on the basic rack and protect with the star screws with a clockwise rotation.

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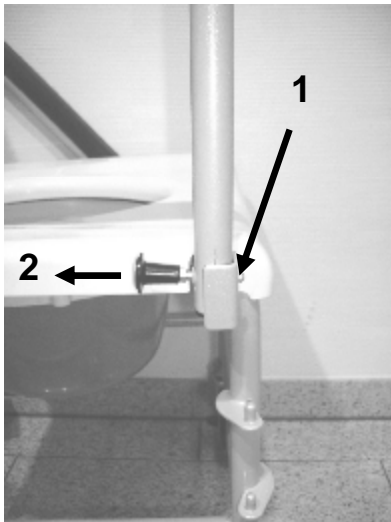
*Pull both armrests upwards until they **audible** engage the interlock.*

The interlock occurs automatically. Make sure that the interlock is successful on both sides, on this occasion, the locking pin (1) must be visible.

In order to unlock the armrests move the locking button to the back (2) and swing the armrest downwards.

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
*To interlock the brake press lever downwards. **This will lock the direction of the wheel in place.***

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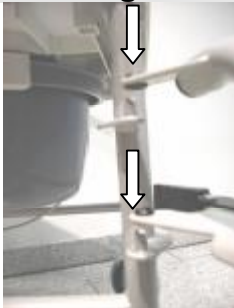
Hang up the provided footrests with both drillings as shown in the respective hold pins

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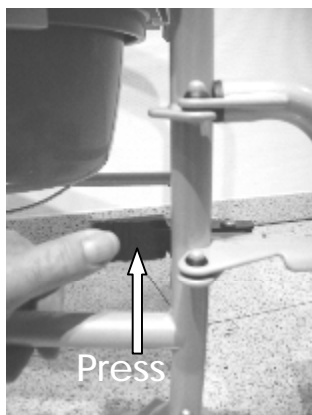In order to interlock the footrests turn them to the middle until its audible engage. In order to unlock press the lever as shown inwards.

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
To adjust the footrests in the height first unlock the fixing clamp (1). After pinching the fixing bolts (2) the footrests can be shifted (3) in 4 steps in the height for leg lengths from 420 – 500 mm.

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ATTENTION: Never use the foot rest as stand up device !

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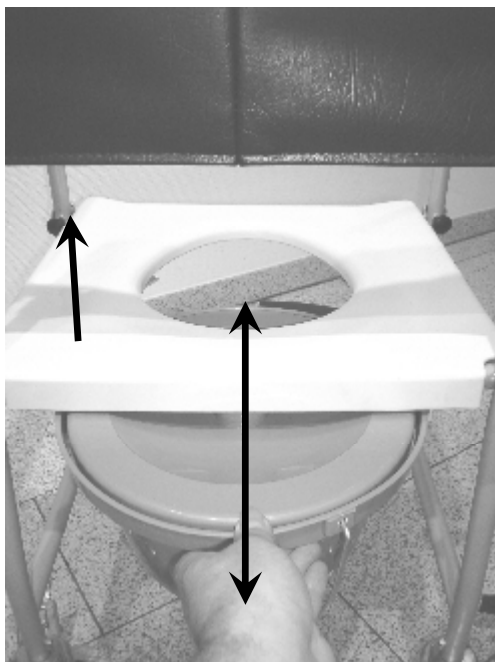Insert or remove the toilet bucket (also with put on lid) always from the back. The toilet bucket engages, the toilet seat can be taken simply upwards for cleaning.

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Clean regularly by hand with a damp cloth and usual household cleaning agents!

Use spraying disinfectants p.e. Descosept AF or Bacillol AF.

Do not use any aggressive detergents!

Screws of TRS 130 must be checked regularly and if necessary retightened.

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Frame: Welded steel tube, chrome plated

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We issue a warranty of 24 months from date of purchase.

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*Subject to technical alterations.
Errors and omissions excepted.*

Stand - as per: **07/2009**

Drive Medical GmbH & Co. KG
Leutkircher Straße 44
D-88316 Isny / Allgäu (Germany)

Drive Medical Limited
Ainleys Industrial Estate
Elland
GB WEST YORKSHIRE HX5 9JP

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